# Kappa light chain Rabbit mAb [O888]

Cat NO. :A93624

## Information:

Applications	Reactivity:	UniProt ID:	MW(kDa)	Host	Isotype	Size
WB,IHC	н	P01601	25 kDa	Rabbit	lgG	50ul,100ul,200ul

### **Applications detail:**

Application	Dilution		
WB	1:1000-2000		
ІНС	1:100,		
ne optimal dilutions should be determined by the end user			

### Conjugate:

UnConjugate

Form:

Liquid

#### sensitivity:

Endogenous

## **Purification**:

Protein A purification

#### Specificity:

Antibody is produced by immunizing animals with a synthetic peptide of Human Kappa light chain.

#### Storage buffer and conditions:

Antibody store in 10 mM PBS, 0.5mg/ml BSA, 50% glycerol (buffer) .

Shipped at 4°C. Store at-20°C or -80°C.

Products are valid for one natural year of receipt. Avoid repeated freeze / thaw cycles.

**Tissue specificity:** 

#### Subcellular location:

Secreted. Cell membrane.

#### **Function**:

V region of the variable domain of immunoglobulin light chains that participates in the antigen recognition (PubMed:24600447). Immunoglobulins, also known as antibodies, are membrane-bound or secreted glycoproteins produced by B lymphocytes. In the recognition phase of humoral immunity, the membrane-bound immunoglobulins serve as receptors which, upon binding of a specific antigen, trigger the clonal expansion and differentiation of B lymphocytes into immunoglobulins-secreting plasma cells. Secreted immunoglobulins mediate the effector phase of humoral immunity, which results in the elimination of bound antigens (PubMed:20176268, PubMed:22158414). The antigen binding site is formed by the variable domain of one heavy

Introduction: WB: Western Blot IP: Immunoprecipitation IHC: Immunohistochemistry ChIP: Chromatin Immunoprecipitation ICC/IF: Immunocytochemistry/ Immunofluorescence F: Flow Cytometry


Cross Reactivity: H: human M: mouse R: rat Hm: hamster Mk: monkey Vir: virus MI: mink C: chicken Dm D. melanogaster X: Xenopus Z: zebrafish B: bovine Dg: dog Pg: pig Hr: horse

For Research Use Only. Not For Use In Diagnostic Procedures.

chain, together with that of its associated light chain. Thus, each immunoglobulin has two antigen binding sites with remarkable affinity for a particular antigen. The variable domains are assembled by a process called V-(D)-J rearrangement and can then be subjected to somatic hypermutations which, after exposure to antigen and selection, allow affinity maturation for a particular antigen (PubMed:20176268, PubMed:17576170)..

# Validation Data:

## Kappa light chain Rabbit mAb [O888] Images

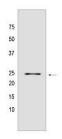

Western blot (SDS PAGE) analysis of extracts from Human tonsil tissue lyaste.using Kappa light chain Rabbit mAb [O888] at dilution of 1:1000 incubated at  $4^{\circ}$ C over

View more information on http://naturebios.com

IMPORTANT: For western blots, incubate membrane with diluted primary antibody in 1% w/v Milk, 1X TBST at 4°C overnight.